### SPECIFICATIONS

additional specifications at zingstanders.com

	Zing Portable	Zing Size 1	Zing Size 2
Weight Limit	36 lbs. (16kg)	70 lbs. (32kg)	154 lbs. (70kg)
Height Range	up to 36" (91cm)	up to 44" (112cm)	40"-60" (102-152cm)
Transfer Height	26" (66cm)	32" (81cm)	31" (79cm)
Hip to Foot Plate Range	8.5"-17.5" (22-44cm)	7"-22" (18-56cm)	15"-33" (38-84cm)
Upper Body Support Height Range	9" (22cm)	4" (10cm)	6" (15cm)
Knee Pad Depth Range	1" (2cm)	5" (13cm)	7" (18cm
Weight of Basic Unit	14 lbs (6kg)	55 lbs (25kg)	82 lbs (37kg)
Frame Footprint	39"x32" (91x81cm)	20"x28" (51x71cm)	26.75"x39.75" (68x101cm)
Degree of Independent Leg Abduction on each leg	15°	30°	30°
Supine (not available on Prone)	0°-60°	0°-90°	0°-90°
Prone (not available on Supine & Supine TT)	0°-60°	0°-20°-MPS, MPS TT; 0°-90°-Prone	0°-20°-MPS, MPS TT; 0°-90°-Prone
HCPCS (MPS, MPS TT)	E0641-Multi-Position Stander	E0641-Multi-Position Stander	E0641-Multi-Position Stander
HCPCS (Supine, Supine TT & Prone)	E0638-Single Position Stander	E0638-Single Position Stander	E0638-Single Position Stander
Positions	Supine	Supine	Supine
	Prone	Prone	Prone
Adjustable Tray	Yes	Yes	Yes
Swing-Away Tray (not available on Prone)	Yes	Yes	Yes
Support	Individual Positioning Components	Individual Positioning Components	Individual Positioning Components
	Padded Strapping	Padded Strapping	Padded Strapping
Dual Control	No	Yes	Yes
Gas Spring Lockout	No	Yes	Yes
Transfer Aid	Near Flat-to-Load	Flat-to-Load	Flat-to-Load
Tool Free Adjustment	Yes	Yes	Yes
Pow'r Up Lift Option	No	No	Yes

#### Hygienic Covers (Size 1 and Size 2 only)

- Available in 5 fun patterns for the pelvic, upper body, lateral, head and tray areas
- Prevents stander from getting soiled, especially during feeding
- · Easily removed and washable











Zing MPS TT-Size 1



ALL SMILES ANIMALS

#### **Frame Colors**



**GREEN** (Portable, Size 1 and Size 2)



PURPLE (Size 1 and Size 2)



(Size 1 and Size 2)




Altimate Medical, Inc. 262 West First Street Morton, MN 56270 USA Tel: 800.342.8968 507.697.6393 Fax: 877.342.8968 507.697.6900 info@zingstanders.com zingstanders.com

Visit our website at zingstanders.com for complete specifications on all Zing products and options.





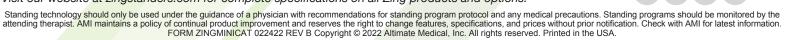

# The Advantage Portable










- The MPS and MPS TT allow supine, vertical and prone positioning without stander
- Smooth transition from a flat-to-load transfer position, supine, vertical and 20 degrees prone - all in one motion, thanks to the unique pivoting action of the
- standers use prone for head control and when fatigue sets in, simply move to

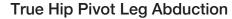

- Pivot points are in line with the hips, promoting proper alignment in the hips, knees and feet
- · Feet maintain the angle with the rest of the legs no wedges needed and no additional adjustment to the knees or foot plates
- Up to 15° (Portable) or 30° (Size 1 and 2) of abduction per leg accomplished through a single knob
- Promotes musculoskeletal development for children during critical growth years
- Abduction can be adjusted and used throughout the entire positioning range

#### Plenty of Growth


- Portable: Fits children from infancy to 36", and up to 36lbs
- Size 1: Fits children from infancy to 44", and up to 70lbs
- · Size 2: Fits children from 40"-60", and up to 154lbs
- No minimum height or weight requirements which makes starting a standing program as young as possible a reality
- · Can easily grow with the child and adapt to their changing needs
- · Many supports, pads and straps are available in multiple sizes
- · Comfortable transfer height with near flat loading makes transfers more manageable as children grow


Zing Portable

#### Remarkably Modular

- · Can be configured as a simple stander or fit a child with the most complex needs
- Designed to grow with the child options can be added as needed

#### **True Multi-Position Stander**

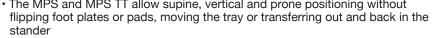

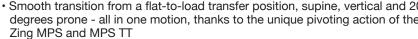

- Provides therapeutic opportunities not possible with other multi-position a slight supine position to rest but remain in the standing position
- · Frame will only go supine and upright with Prone Position Lockout engaged



Supine Lying - Flat-to-Load

- The MPS, Supine, MPS TT and Supine TT feature a multi-use tray that is angle adjustable and swings away either direction
- Tray can remain attached during transfers
- More activities are possible by swinging tray away during standing therapy

#### **Unique Lift Mechanisms**

- Dual controls are available allowing angle adjustment by foot pedal or hand control
- · Gas Spring Lockout prevents an inadvertent change of position, especially important when around other children
- The Pow'r Up Lift option on Size 2 brings the user to the standing position with the touch of a button

#### **Accommodating Contractures**

- No special modifications are needed to accommodate children with contractures
- Using the multi-adjustable foot plates, multi-adjustable knee pads and the independent calf pads, it is possible to adjust for contractures in one or both sides

#### Easy Adjustment

- The Portable offers built in straps on the support pads
- · Tool-free adjustment of all components
- · Every adjustment can be made with a simple turn of a knob or ratchet handle

#### **Full Line of Standers**

- Zing Portable: MPS Supine Prone Zing Size 1: • MPS Supine Prone
  - MPS TT Supine TT
- Zing Size 2: Supine Prone

MPS TT
Supine TT

## SUPINE TO PRONE







Zing MPS TT-Size 2



Zing Portable

Zing MPS-Size 1

Great for a multi-user facility or school

